



Jan. 2020

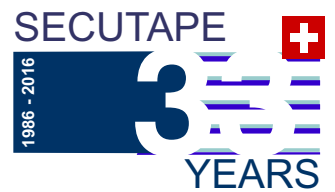

FIXATIONS

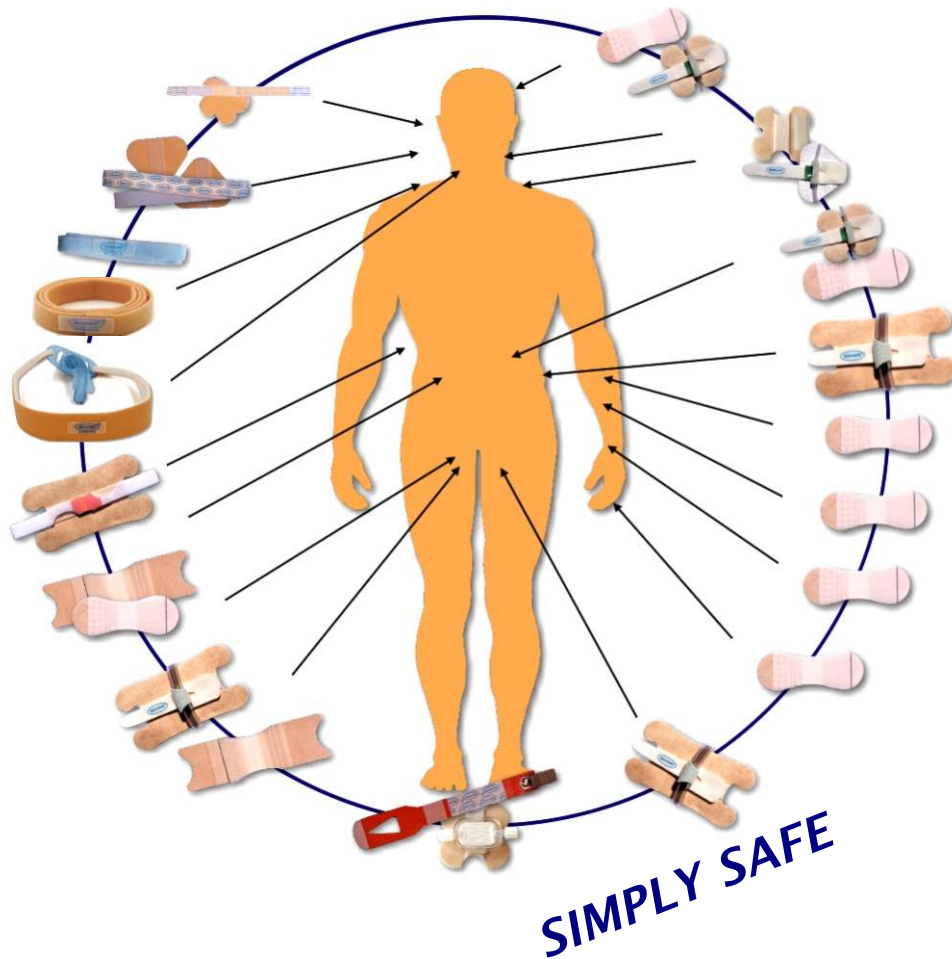

CATALOG



TechniMed AG

SWITZERLAND since 1987
fon 0041 / 71 844 49 20

Kirchstrasse 54
CH-9400 Rorschach

fax 0041 / 71 845 52 45

Germany

fon 0800 007 9892

fax 0800 941 3012

mail: info@technimed.ch

Austria

fon 01/ 440 7979

fax 01/ 440 7990

CE

Made in Germany

www.secutape.ch

Copyright by TechniMed AG


The SECUTAPE idea:

The base plaster adheres completely to the skin; there are no detached spots!

The 'tube' is embedded in the adhesive and a velcro prevents the Flap from opening!

Via this Flap the movement of the tube is transferred to the whole plaster.
The fixation holds and can be opened and closed again, using the velcro flap.

All SECUTAPE fixations working in the same way which we have patented 1986.

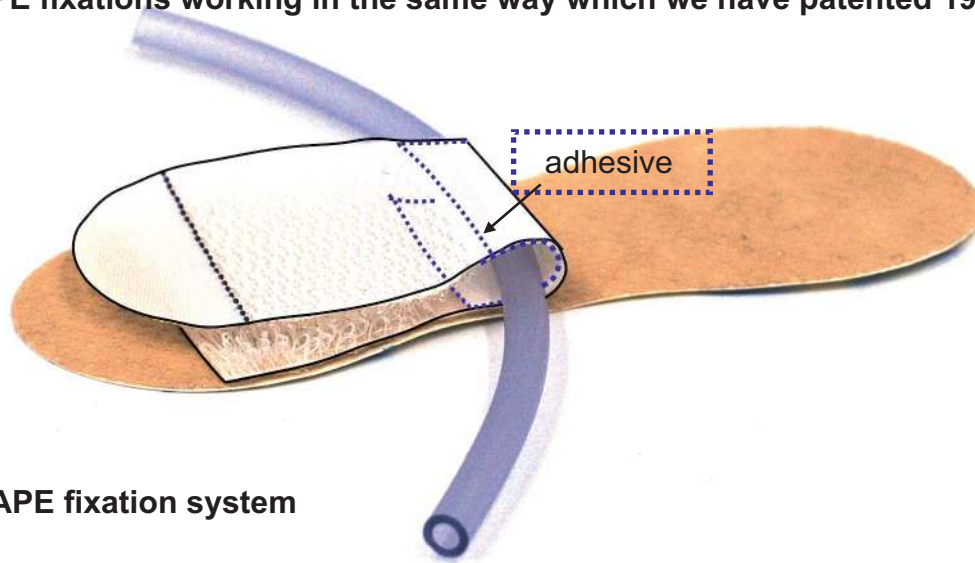

The SECUTAPE fixation system

Provides safe and reliable solutions for practically all attachable catheters, drains, probes, tubes and cables; i.e. for all tube-shaped accessories.

We are the specialists in this field. For thirty years we have been manufacturing and distributing exclusively millions of fixations and we hold quite a number of patents.

We are offering our customers free training and advice through our field workers, including extensive documentation.

Our INFOCARE E-Learning platform for training of new employees and to refresh existing knowledge.

Our fixations are high-tech products that are manufactured in Germany. The materials used are manufactured specially for us and come from Switzerland, Germany, the USA and the UK. They are hypo-allergic, medically safe and contain no solvents.

Only through the use of high-quality materials in combination with our special and ergonomically designed shapes, we have achieved the astonishing holding power of our fixations.

Our products certainly come at a cost. But they hold better, do not cause skin damages and must be changed less frequently than others.

Calculated over the period of use, they are not expensive after all.



Originalgrösse



SECUTAPE®

Size 2 / 2+ / 3+

Standard fixation


Since 30 years our reliable fixation!

for the fixation of :

- PEG
- CAPD
- CVC
- suprapubic catheters
- infusions
- urin-catheter
- redondrains
- ports
- wounddrains
- and all other tubes

For heavy tubes pls. use size 3+

The fixation can be opened and closed many times.

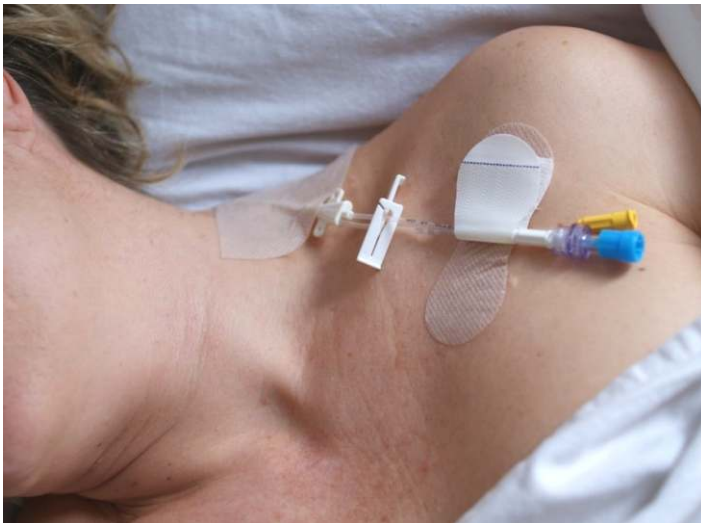

ord.no. 8500200 1 PU 200 pcs.
foil basic
singlepack 8 bags a.25pcs.

ord.no. 8500210 1 PU 1000 pcs.
20 bags. à 50 pcs.



ord.no. 8500201 1 PU 200 pcs.
NW basic
singlepack 8 bags a.25pcs.

ord.no. 8500211 1 PU 1000 pcs.
20 bags. à 50 pcs.



ord.no. 8500301 1 PU 200 pcs.
NW basic
singlepack 8 bags a.25pcs.

ord.no. 8500311 1 PU 1000 pcs.
20 bags. à 50 pcs.



Stays appr. 3 - 5 days

Info

ord.nb. 8500200 sticks better to bandaid,

ord.nb. 8500201 stays better on the skin.



Made in Germany
www.secutape.ch



Originalgrösse



SECUTAPE®

Velcro-binder

**A new universal system
for all applications and diameters!**

**Velcro binder holds reliably, without any adhesive
tubes and wires with the new patented NANOPLAST®
antislip foil.**


The fixation can be **opened and closed many times**
without leaving any marks on the tube!

Velcro binder comes as a set with different holders,
also in Hydrocolloid for every application.

Due to the high quality breathing material, the
fixation stays even when it gets wet.

NANOPLAST® has the so-called Gecko effect
and is made out of silicon with aprox. 29000 nano legs per cm .


**The fixation is very good for urine catheter
and especially for the fixation of implanted tubes/lines.**

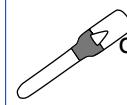

Velcro-binder
ord.no. 8501000 PU 200 pcs.
singlepack. 8 bags a 25 pcs.





Velcro binder - Set SMALL
Hydrocolloid
ord.no. 85001100 PU 100 pcs.
singlepack. 4 bags a 25 pcs.

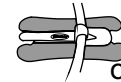

Velcro binder - Set BIG
ord.no. 8501130 PU 100 pcs.
singlepack 4 bags a 25 pcs.

**Stays up to 3 days,
Hydrocolloid up to 5 days**

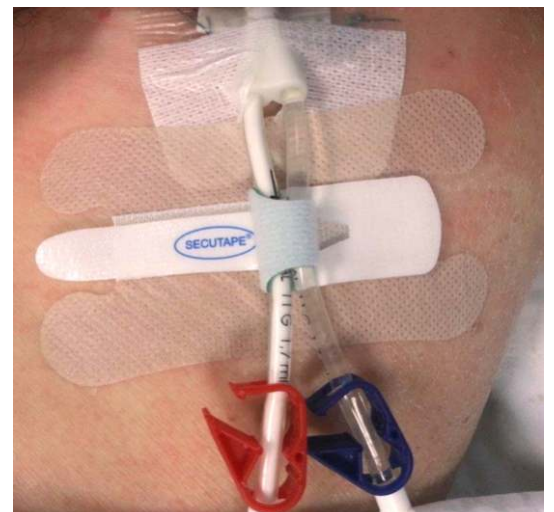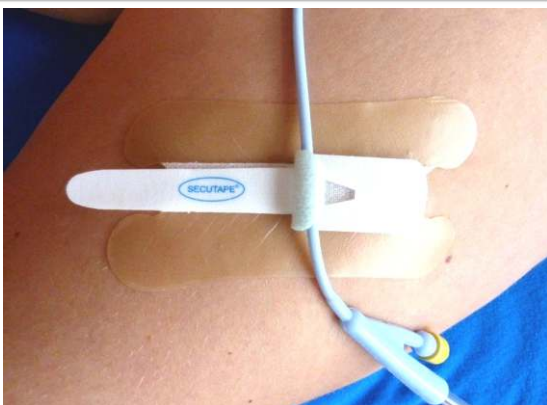


CE

Made in Germany

www.secutable.ch



Originalgrösse

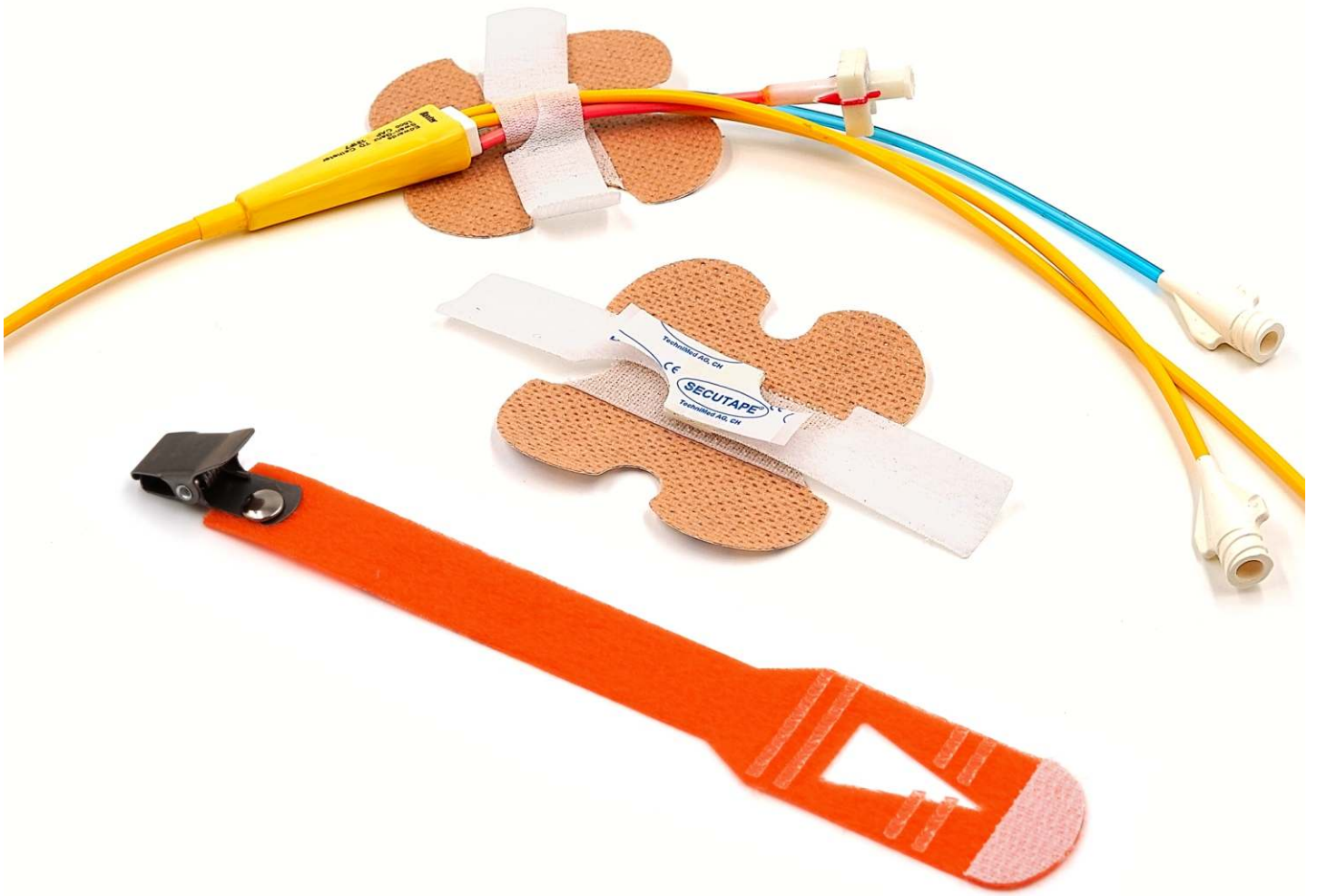



Fixation for "tubes and wires" especially for CVC

UNIFIX is perfectly suited for the fixation of thinner catheters, Probes and lines, e.g. multi-lumen CVC...

The basis is a very good holding, ergonomic shaped holder made of skin-friendly, breathing tape that also holds when it gets wet.


The catheter is placed in an adhesive bed and held securely with a Velcro strap.




Fixing Set CVC

The set includes a UNIFIX and a non-slip flexible holder that can be attached to clothing or bed linen, as an **additional fixation**.

The use of two fixations reduces strain and movement on the CVC, which prevents infections.




Unifix
ord.no. 8501500 PU 200 pcs.
8 bags. à 25 pcs. singlepack




Fixing set CVC
ord. no.. 8501510 PU 100 pcs.
4 bags. à 25 pcs. singlepack

UNIFIX can also be used to attach filter, transducer or similar parts securely to the body.

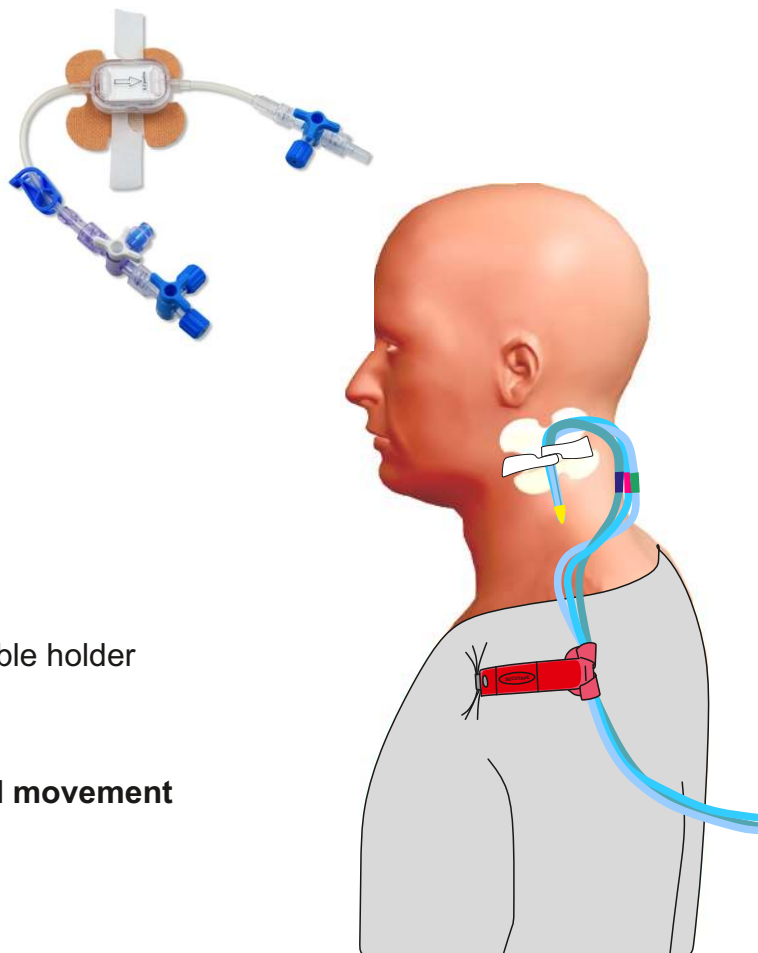



Originalgrösse

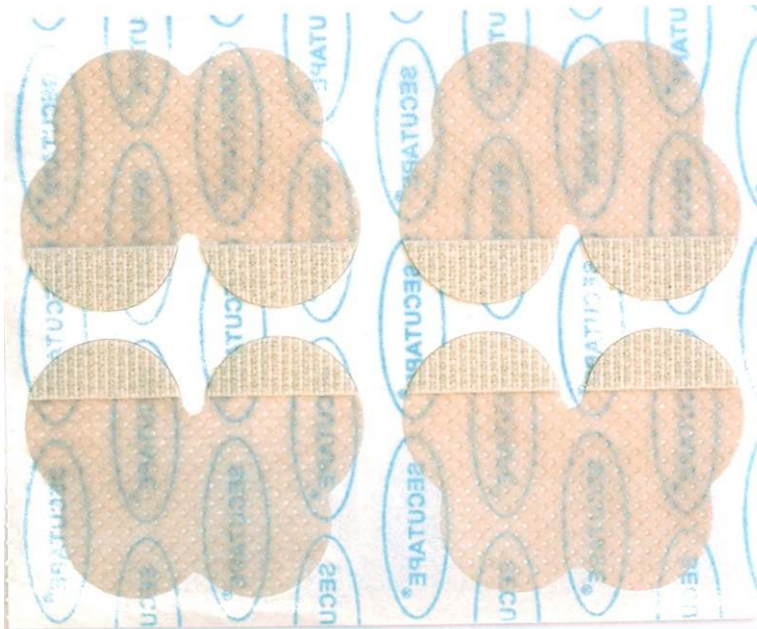



Nasal fixations

Fixation for nasale sondes

Patented two piece fixation for all nasal sondes.

A specially shaped basic holder made of a skin friendly and breathable nonwoven or a hydrocolloid.

The basic holder has a Velcro part, a second strap with Velcro is glued to the tube.

The holder remains on the nose, if necessary only the straps are changed.

There is no more skin damage caused by the constant tearing of the plaster strips. The probe is hanging loose in the nose without causing any problem.




ord.no.. 8500400 1 PU 200 pcs.
hydrocolloid



singlepack.
set




ord.no.. 8500401 1 PU 200 pcs.
nonwoven



8 bags a.25pcs.
singlepack.
set




ord.no.. 8500410 1 PU 400 pcs.



board à 4 pcs.
singlepack.



ord.no.. 8500440 1 PU 500 pcs.




board à 5 pcs.
singlepack.

Stays appr. 3 days
hydrocolloid 5 days



Originalgrösse

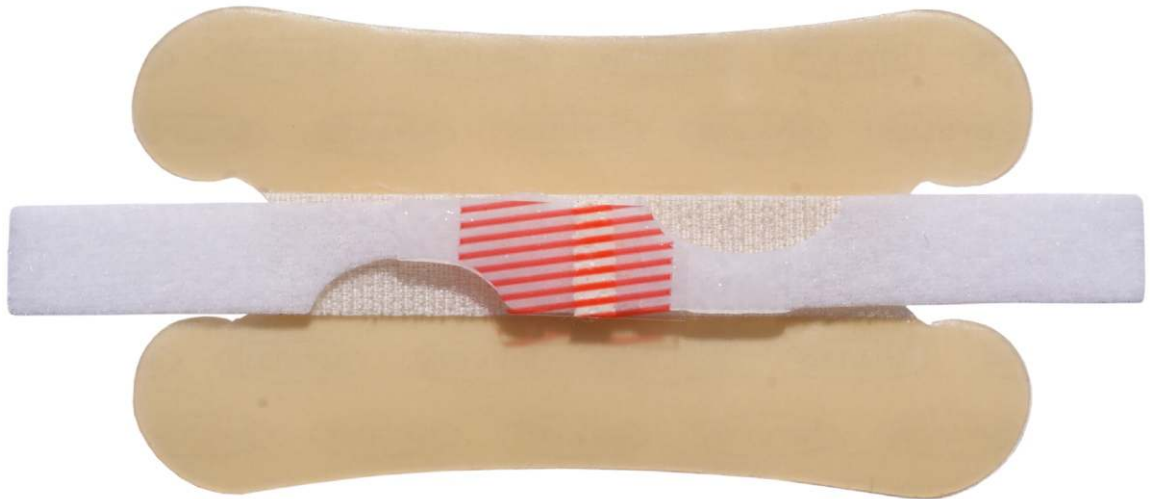



Fixing set for big lumina

Our best holding Fixation for large tubes!

Good for mobile patients!

Specially designed for:
chest tubes, ECMO, rinsing catheters and
for the fixations of urinary catheters.

Through a special coating it even holds
silicon catheters!

Due to the high quality nonwoven material,
the fixation stays even when it gets wet.



ord.no. 8501200 1 PU 100 pcs.
hydrocolloid 4 bags. à 25 pcs.
singlepack.



ord.no. 8501201 1 PU 100 pcs.
nonwoven 4 bags. à 25 pcs.
singlepack.



**Stays appr. up to 4 days
hydrocolloid up to 6 days.**



Originalgrösse



SECUTAPE®

Fixing set for endotracheal tubes oral/nasal

Since over 30 years our reliable fixation for endotracheal tubes!

Easy to fit by just one person and inexpensive due to longer application time..

Two holders of a very good adhering and skin friendly nonwoven or hydrocolloid with Velcro

Non stretchable, hypoallergenic nylon straps which do not absorb liquid.
The straps are centrally coated with a special adhesive.

The system holds very well when used correctly and allows safe aspiration without a second person.

The standard strap is **soft on one side**.

The Softband is **soft on both sides**, which prevents skin damages.



**Stays appr. up to 3 days,
hydrocolloid up to 5 days.**



Originalgrösse

Band 40 cm



CE

Made in Germany

www.secutape.ch

SECUTAPE®

Tube fixing set EXTRA


Our known tube fixing set with a newly developed strap.

With our proven nonwoven holders and a very soft hypoallergenic nylon strap that does not absorb humidity and does not stretch.

The band has a slot in the center for a secure position of the tube.


The band is very well held on the tube with a self-adhesive Velcro strip.

The system holds very well when correctly installed and allows safe aspiration without a second person.

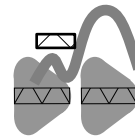

Fixing set EXTRA for endotracheal tubes

ord.no. 8500610 1 PU 200 pcs.
singlepack 8 bags à 25 pcs.
Both parts are available separately

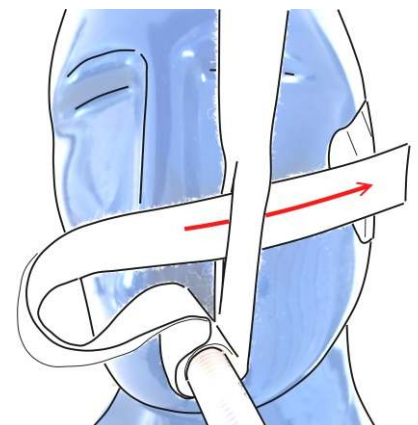

CE

Made in Germany
www.secutable.ch



Länge 37 cm



Liegedauer je nach Anwendung:
bis 3Tage
Hydrocolloid - 6 Tage

CE
Made in Germany
www.secutape.ch



Short-term fixing tape

self-adhesive
for endotrachealtubes.

Fixing tape for the short-term fixation of endotracheal tubes.

Newly developed high-quality self-adhesive fixing tape,
specially designed for the anesthesia.

The hypoallergenic, skin-friendly and excellent
holding tape does not stretch and rip.

It is packed in a hygienic dispenser and easy to handle
due to its shape and the special liner.

The perforated dispenser is opened at the front, which
eases the withdrawal and storage.



Short-term fixing tape

for endotracheal tubes

ord.no. 8500820 1 PU 400 pcs.
box a. 100 pcs.

ord.no. 8500821 1 PU 2000 pcs.
box a. 100 pcs.




**This strap does not replace
a good long-term fixation!**



CE

Made in Germany
www.secutape.ch

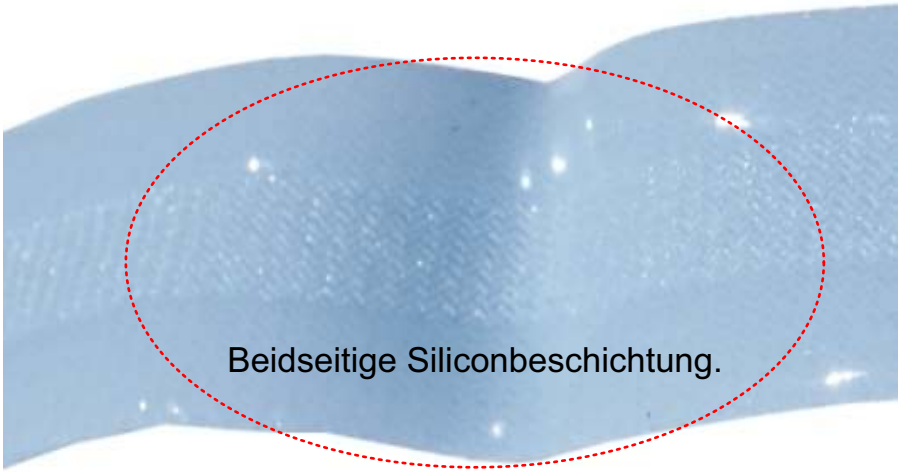

Beidseitige Siliconbeschichtung.

SECUTAPE®

Tube-band


Pat. fixation for endotracheal tubes

A hypoallergenic soft nylon band, it does not stretch and does not absorb any liquid.

Since the strap does not stretch no need to tight.

Due to the silicon coating on one side the strap has no hard edges, does not slip on the skin and **the tube stays in place!**

It is easy in handling, prevents most skin problems and must be changed less frequently than conventional textile or foam straps.

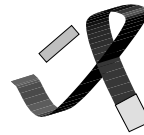

Tube-band

ord.no. 8500850 1 PU 200 pcs.
singlepack. 8 bags à 25 pcs.
75 cm long



Stays approx. up to 2 days.



CE

Made in Germany
www.secutape.ch



Länge 75 x 2cm





Universal-band

For endotracheal tubes & tracheal cannulas

Newly developed double-sided nylon band is hypoallergenic and pleasantly soft.

It does not stretch and does not absorb liquid. Therefore the patient does not have a humid band around the neck.

The band is also easy to clean.

It has to be changed less frequently than the usual foam or fabric bands.

Since the band should not be applied tight, there will be less skin damage and more comfort for the patient.

One product for both applications!



Time of use approx. up to 3 days.





Originalgrösse

Tracheofix
53 x 3.5 cm

Superflausch



Superabsorberpad

Masse 21 x 4.5 cm





Tracheofix

Fixation for tracheal cannulas

Newly developed thick and fluffy nylon band.
Hypoallergenic and pleasantly soft.

The high quality band does not stretch and does not absorb liquid.
It is easy to apply and has to be changed less frequently than the usual foam or fabric bands.

For long term fixation or patients with skin problems.
Use the additional super absorber pad.

The medical grade absorber keeps the skin dry and prevents skin damage.

Pull the liner off the pad and stick it inside the band.



Tracheofix
LARGE length 53 cm
 ord.no. 8501700 1 PU 100 pcs.
 singlepack 8 bags à 25 pcs.

SMALL length 40 cm
 ord.no. 8501720 1 PU 100 pcs.
 singlepack 8 bags à 25 pcs.



Superabsorber Pad
 length 20 cm
 ord.no. 8501770 1 PU 100 pcs.
 singlepack 8 bags à 25 pcs.



Stays approx. up to 5 days





Länge 50 x 3 cm
Absorber 20 cm



CE
Made in Germany
www.secutape.ch



Fixation band for tracheal cannulas


A thick and fluffy two-piece Nylon band.

The high quality hydrophob and hypoallergenic band is easy to apply, it does not stretch and does not absorb liquid.

Due to its material, it needs to be changed less frequently than conventional foam or textile bands.

It fits all sizes, just cut it off if it is too long.

Stays up to 3 days

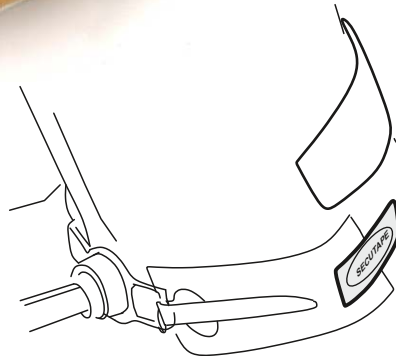

Tracheofix Comfort Set


The band is also available as a set together with a medical grade super absorber pad.

It is a good solution for patients with skin problems or very sensitive skin.

The enclosed med. grade absorber pad keeps the skin dry and is a protection against skin damage.

Therefore it is an ideal long-term fixation and particularly suitable for home care.

Stays up to 5 days.

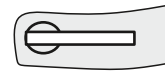

Tracheofix Comfort
length 50cm

ord.no. 8501750 PU 100 pcs.

singlepack. 4 bags. a 25 pcs.



Tracheofix Comfort-Set
With super absorber



length 45cm, absorber 20 cm

ord.no. 8501740 PU 100 pcs.

singlepack. 4 bags. a 25 set



Made in Germany
www.secutape.ch

SECUTAPE®

Superflausch



PVC & LATEX FREI!

Fixierhalter zum Einsatz bei:

Tubing

Hand Art. 8500800

Art. 8501000

Reinigen des Halters
nur mit Wasser und Seife reinigen!

Nicht sofort belasten!!

Hydrocolloid: braucht Hauttemperatur um auf zu halten!

Weitere Informationen für
in unseren Anwendungs-
oder unter
www.secutape.com.

CE

Made in Germany
www.secutape.ch

SECUTAPE®

Head holder-set ECMO

Adjustable head holder set for fixation of an veno-arterial ECMO

The elastic head holder is a thick and soft nylon band and fits all head sizes.

The holder is, in combination with the Velcro straps and the hydrocolloid holders, a very reliable Fixation, which gives the patient a certain room for movement.

The additional Fixation on the neck helps to reduce the strain on the hoses.

With the enclosed KROKOFIX, which should be fixed to the bed linen, you reduce the strain on the drains.



CE

Made in Germany
www.secutable.ch



Originalgrösse






Krokofix

Fixation for "tubes and wires" at clothes and bed linen.

- two functions!
- simple handling !
- two sizes!
- secure and non slipping adjustment!
- variable use!
- disinfection with usual means!

A manual anti slip fixation for tubes of all sizes.
Good for mobile patients and strong pulling.


Just clamp it to clothes or bed linen.

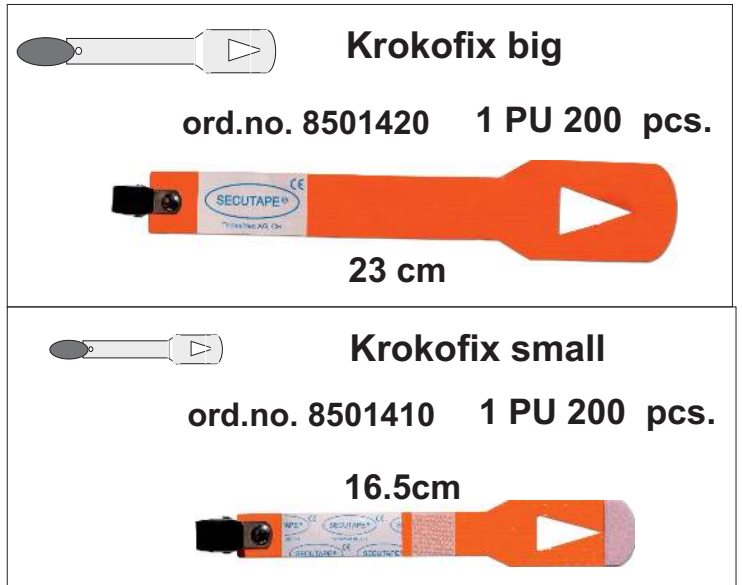

SECUTAPE®

**FIXIERUNG
FIXATION
FISSAGGIO**

SIMPLY SAFE

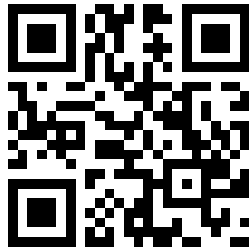

www.secutable.ch